



Research on Surface Sampling for Determination of Pesticide Residues in Pome Fruit

Yunxia Luan, Hua Ping, Zhihong Ma, Ligang Pan

► **To cite this version:**

Yunxia Luan, Hua Ping, Zhihong Ma, Ligang Pan. Research on Surface Sampling for Determination of Pesticide Residues in Pome Fruit. Daoliang Li; Yingyi Chen. 5th Computer and Computing Technologies in Agriculture (CCTA), Oct 2011, Beijing, China. Springer, IFIP Advances in Information and Communication Technology, AICT-369, pp.516-527, 2012, Computer and Computing Technologies in Agriculture V. <10.1007/978-3-642-27278-3_53>. <hal-01361025>

HAL Id: hal-01361025

<https://hal.inria.fr/hal-01361025>

Submitted on 6 Sep 2016

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Distributed under a Creative Commons Attribution 4.0 International License

Research on surface sampling for determination of pesticide residues in pome fruit

Yunxia Luan^{1,2}, Hua Ping^{1,2}, Zhihong Ma^{1,2}, Ligang Pan^{1,2}

¹Beijing Research Center for Agrifood Testing and Farmland Monitoring, Beijing, China

²National engineering research center for information technology in agriculture, Beijing, China

{ Yunxia Luan, luanyx@nercita.org.cn

Abstract. A surface sampling method for analysis of pesticide residues in pome fruit was studied by layer-by-layer scan from pericarp to kernel. Phoxim residues of apples, pears and papayas was determined by GC-MS, respectively, in core, different layers of pulp, peel, as well as the whole fruit. Phoxim residue content was decreased when depth of the layer increased. It was inferred from experimental data the the appropriate depth of surface sampling is 2.5mm, 2mm and 3.5mm in apples, pears and papayas. Although the ratio of the pesticide in pericarp and the total (Cp/Ct) ranged from 6.241 to 9.262 in apple, 4.254 to 5.069 in pears, and 4.39 to 5.037 in papayas in different groups, the ratio of Cp/Ct and quality of specific group is essentially unchanged for a certain kind of fruit. According to the results, a conversion formula of the content of pesticide residues in whole fruit and pericarp was set, based on surface sampling method.

Key words: Poem fruit, Pesticide residues, Phoxim, Surface sampling

1. Introduction

One of the major roles of public health agencies is to ensure safe products for consumers through analysis of residual agricultural chemicals and veterinary drugs in foods. The use of agricultural chemicals and veterinary drugs in agriculture is necessary to improve the quality of the food produced. They play a beneficial role in providing a plentiful, low-cost supply of high-quality food. On the other hand, as a consequence of this use, the presence of residues in food that is a critical element of overall public health is unavoidable, and residues in food are of great importance in the evaluation of food quality. Organophosphorus (OP) pesticides are the most widely used agricultural pesticide. More than 70% of that are hyper-toxic or high-toxic, and normally are forbidden to use in planting fruits and vegetables. OP insecticides are routinely applied to fruit crops (e.g., apple, peach, pear) that are subject to consumption as single servings [1]. Consumption of individual fruits containing high OP insecticide residues may result in high exposure over a short period of time. Residue testing for the presence of pesticides, however, is routinely performed on composites of 5 to 10 individual or unit samples rather than samples of individual fruit [2].

Experimental study of the fruit for pesticide residues in agricultural products exist in the form, that is attached to the fruit into the fruit surface and the distribution of pesticides within the fruit. It is necessary to provide a representative sample to the laboratory for analysis [3]. No matter how carefully

laboratory analyses are carried out, the results will be of limited value unless the sample provided is truly representative of the sample being tested [4]. The uneven distribution of pesticide residues among the treated objects leads to an inevitable variability of pesticide residue levels measured in the samples, which may significantly contribute to the combined uncertainty of the analytical results. Data generated in the UK have indicated that pesticide residue levels can be highly variable between the individual fruit from the same crop or lot in trade. Many fruit or vegetables in trade are mixed after harvest to form combined lots [5]. Analysis of composite samples provides a good indication of average residues but, where the lot has been mixed, such average values are potentially misleading [6]. Residues monitoring is the best means available for general control of pesticide use and consumer exposure, but new strategies for sampling and analysis are required to address the combined effects of residues variability and mixing of lots.

Development of analytical methods for pesticides in agri-food requires precise method of sampling. Pome fruit peel is the main remaining parts, parts of the peel residue is a measure of its ability to meet the safety standard limit fruit key factor. The apple pesticide residue question has become the chief obstacle in affecting the apple industry internationalization of our country [7]. At present, samples of pome fruit pesticide residue mainly from the entire apple or some edible parts, because the results from these have the generalization. However, the pesticide residue mainly exists in the pericarp. Therefore, the careful research to pesticide residue of the pericarp provides the theory instruction and the design basis to the pesticide residue and the elimination research.

However, limited data have been reported that directly compare the sampling procedures against one another. Even less information is published regarding the most appropriate part of the material to use in conjunction with a sampling procedure for pesticides detection, which is particularly important to the portable rapid detector of contamination in fruits and vegetables [8]. In our Previous studies, the temporal and spatial variations of phoxim residue in apples were studied. The results show that half-life of phoxim in apples was 1.64 days. The phoxim residue order in different position of apple is: pericarp >entire apple> sarcocarp> kernel and the ratio of residue in pericarp and the entire apple ranged from 6.241 to 9.262 in groups with different specific surface area [9]. The results indicate that, there are significant differences in the contents of residues between the pericarp and entire apple, which provide evidence in the theory instruction and design basis to the sampling method of pesticide residue based on pericarp. In this paper, further research was carried out on a surface sampling method for analysis of pesticide residues in pome fruit was studied by layer-by-layer scan from pericarp to kernel, and to explore the relationship between the content of pesticide in pericarp and the whole fruit.

2. Materials and methods

2.1. Chemicals and Reagents

Phoxim standard was purchased from National Institute of Metrology P.R.China. High-purity solvents cyclohexane, acetone, and acetonitrile were all purchased from Fisher Chemicals [10]. Analytical reagents magnesium sulfate and sodium chloride were purchased from Beijing Chemical Reagents Company (China). High-purity solvents acetonitrile, suitable for liquid chromatography, gas chromatography, and residue analysis used in the extraction and cleanup steps, were purchased from EMD Biosciences Inc(Mississauga, ON). Reagent grade $MgSO_4$ and sodium chloride were also purchased from EMD Biosciences Inc.

2.2. Apparatus

The gel permeation chromatography (GPC) equipped with gas chromatography mass spectrometry system (GC/MS) was used. The GPC consists of two LC-10ADvp pumps, a SIL-10ADvp auto-sampler, a Shodex CLNpak EV-200AC column (2mm i.d.×150 mm) and CTO-10ASvp column oven, a SPD-10Avp UV detector, two FCV-12AH flow channel selection valves (RV.A, RV.B) and a SCL 10Avp system controller. GC/MS machine is a Shimadzu GC/MS-QP2010 instrumentation equipped with a PTV-2010 large-volume injection device.

GC/MS data analysis was triggered by a contact closure start signal from the HPLC controller. Data acquisition was performed using a C-R8A plus data processor. All these parts are the products of Shimadzu (Kyoto, Japan), except the Shodex CLNpak EV-200AC column (Shoko Co., Tokyo, Japan). Acetone/ cyclohexane mixing solvent (3/7, v/v) was used as the mobile phase of GPC [11], and the flow rate was set at 0.1 mL/min. The mobile phase was degassed using DGU-14A degasser (Shimadzu), and the GPC column was kept at 40 °C in the column oven.

2.3. Sample preparation

All samples were purchased at local markets in Beijing. Samples were extensively collected to achieve good sample homogeneity. The phoxim compounds were applied using an air blast sprayer and application rates guaranteed to contain 50% phoxim emulsifiable with the dilution multiple of 1:500. Individual fruits were cored and sliced into 5 segments (pericarp(P), sarcocarp1(S1), S2,S3,S4 and kernel) using a corer/slicer retailed and peeling machine for domestic use. The first and alternate slices of each apple were taken, chopped manually using a knife, and placed in a plastic bag for storage at -25 °C until extraction and analysis. The remaining segments of individual apples were retained in separate bags for preparation of composite samples. Composites were constructed from the apples prepared from each group by randomly selecting 8 apples for each treat, without regard to the style and weight of each apple. Composite samples were prepared by thoroughly mixing the chopped apple pieces from the bags containing the retained portions of the individual apples. Composite samples were frozen until extraction and analysis. To prepare each sample, 10 g of a previously homogenized food material was transferred into a suitable glass vessel. Then, 10mL acetonitrile was added to each sample using an adjustable-volume solvent dispenser. The glass vessels were capped before vortex mixing for 1 min at maximum speed. Once the initial sample mixing was completed, 1 g NaCl and 4 g anhydrous MgSO₄ were added and mixed immediately on a Vortex mixer for 1 min. It is important to note that this step must be taken immediately after the initial mixing step to prevent the formation of MgSO₄ conglomerates. To separate phases, samples were centrifuged for 10 min at 1570×g. Using an adjustable repeating pipette, 1.0mL aliquot of upper acetonitrile layer was transferred into a 1.5mL flip-top microcentrifuge vial containing 150 mg anhydrous MgSO₄ and 50 mg PSA sorbent. The vial was tightly capped and shaken on a vortex mixer for 1 min before extraction. Then the mixed extraction solution was centrifuged for 5 min to separate solids from solution. The solution was then transferred into an autosampler for GPC-GC/MS analysis.

2.4. Determination of phoxim in apple by GPC-GC/MS

Analysis was performed by GPC-GC/MS. The analytical GC/MS equipped with a PTV-2010

large-volume injection device was carried out using deactivated silica tubing (5m×0.53mm i.d.) , a Rtx-5MS pre-column (5m×0.25mm) and a Rtx-5MS column (30m×0.25mm×0.25µm film thickness; coated with 5% phenyl and 95% methylpolysiloxane; Restek Corporation, Bellefonte, Palo Alto, USA). The temperature of the PTV injector was set at 120°C for the initial 5 min of sampling time, and then increased to 250°C at 100°C/min. The oven temperature was maintained at 80°C for 5 min, subsequently increased to 280°C at a rate of 8°C/min, and then held constant for 10 min [12]. The quadrupole mass spectrometer was operated in the electron impact ion (EI) mode. Ion source temperature and interface temperature were set at 200°C and 250°C, respectively. The mass spectrometer was operated in an ionizing energy of 70 eV. Injection volumes were 1 µL for all analyses. Helium (99.999%) was used as the carrier gas [13].

2.5. Quality Assurance/Quality Control

Apples from market with no phoxim insecticide application were obtained, prepared the same as the rinsed and peeled apples, and were used to prepare blank apple matrix for quality assurance testing. With each set of samples 10-15 samples per set) extracted and analyzed, two aliquots of blank apple matrix were prepared for separate extraction followed by extraction, cleanup, and analysis as for all other samples. Background levels of some analysis were periodically detected in the blank matrix samples and were used for background subtraction in the determination of recovery from spiked matrix only; residue concentrations in samples were not blank corrected. No traces of the phoxim were observed in any of the reagent blanks analyzed.

3. Results and discussion

3.1 Determine of the parameters for surface sampling

In order to determine the exact thickness of pericarp and sampling style, the material was cut by 1.5mm-2mm each piece from the pericarp to kernel, and pesticide residues in fruit peel, kernel and pulp of different depth were determined. Taking the Fuji apple as example, according to dilution ratio of 1:1000 dilution of Phoxim EC, more than fifteen days after spraying pesticide on fruit, pesticide residues in samples was detected according to the method of NY/T 761.1-2004 "vegetables and fruits organophosphorus, organochlorine, pyrethroid and carbamate pesticide multi-residue methods-Part 1 of fruits and vegetables organophosphorus pesticide multi-residue method" and GB/T5009.199-2003 "vegetables of organophosphorus and carbamate pesticide residues in rapid detection." for lab and rapid detection. Each sub-sample was taken from 15 samples of different shape and weight. The results are presented in Figure 1.

A

B

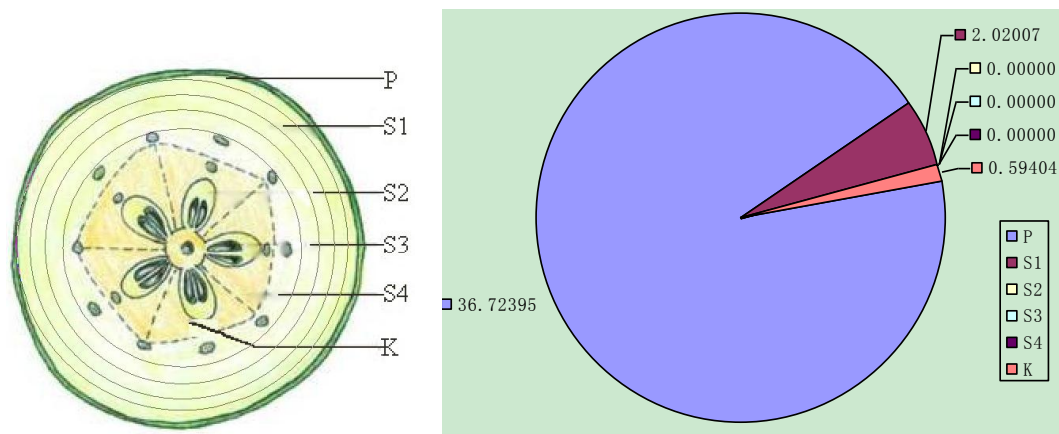


Figure1. Compartmentalization and sampling method of apple (A) and the content of Phoxim in each depth (B)

P: pericarp of apples with 0.5mm pulp, S1: The first layer with 2mm pulp, S2: the second layer with 2mm sarcocarp, S3: the third layer with 2mm sarcocarp, S4: the sarcocarp near kernel, K: kernel or core of apples.

Experimental study of the pesticide residues in pome fruit, which is attached to the fruit surface and the distribution of pesticides within the fruit was carried out. From experiments performed and analysis of pesticide residues in different layer from pericarp to kernel, it is showed that, pesticide was detected mainly in layer P and S1, the pesticide in pericarp and the first layer with 2mm pulp account for 98.49% of the whole sample. Through different parts of the fruit residue test results revealed that the major part of remaining in the apple peel, fruit pulp and other parts of the residual is minimal. According to the apple skin to the stone from the stratified analysis, pesticide residues, mainly in the skin and pulp of the following 2.5mm, this part of the sample selected sampling to detect the levels of residues of pesticides, because this part of the high levels of pesticide residues, can increase rapidly screening test accuracy.

According to the experimental results, the standard design of the skin sampling methods, reduction in access to sub-samples the following 2.5mm thick skin and flesh as the test object.

3.2 Study of the Weight Coefficient of pesticide residues in pericarp

Selected based on the standard higher levels of pesticide residues in samples of skin as a target for different types of fruits, pesticide residues on its surface and overall there is a certain amount of pesticide residue conversion relationship [14]. The research was supported by the skin of pesticide content in the trend over time to determine the appropriate sampling time interval after application, and in accordance with the standard parameters of the specific sampling techniques to obtain test samples, the research team through the following experiments, made pome fruit surface pesticide residue content and the overall coefficient of pesticide content in the skin pesticide residues.

Using Apple, pear, and papaya(n=15), sampling the surface homogenized pericarp(with 2.5mm pulp) as detection material, homogenized with a whole apple parallel processing [15], determination of pesticide content, access to pesticide skin factor, diluted with water according to the proportion of 1:1000 dilution phoxim EC, using the method of spraying pesticide on apples [16]. According to

NY/T761.1-2004 "vegetables and fruits organophosphorus, organochlorine, pyrethroid and carbamate pesticide multi-residue methods-Part1 of fruits and vegetables organophosphorus pesticide multi-residue method" approach with gas chromatography of whole fruit and pericarp of the accurate detection of pesticide residues pesticide residues in order to obtain skin factor [17].

As the result in Fig2-4, the three kinds of fruit showed various levels of pesticide in pericarp and the whole fruit, though the content in pericarp was much higher than the total homogenate. So, the ratio of the pesticide in pericarp and the whole fruit (Cp/Ct) was got from the three kinds of fruits with different shape and quality as shown in Table1.

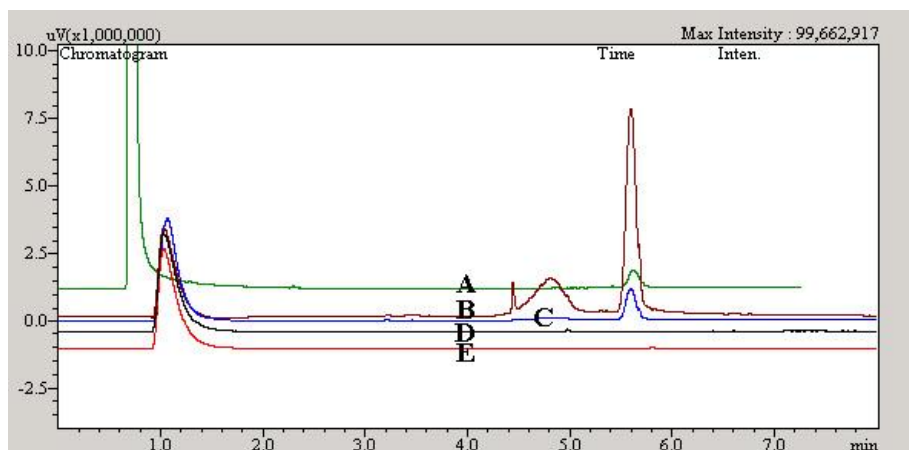


Figure 2: GC diagram of phoxim in Fuji apple

A phoxim 5.25 mg/kg; B sample of pericarp; C sample of the whole fruit; D and E as CK

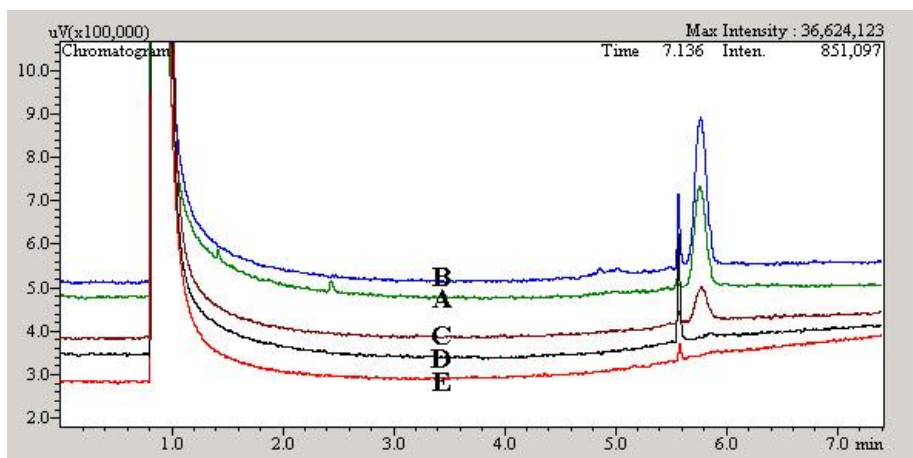


Figure 3: GC diagram of phoxim in pear

A phoxim 5.25 mg/kg; B sample of pericarp; C sample of the whole fruit; D and E as CK

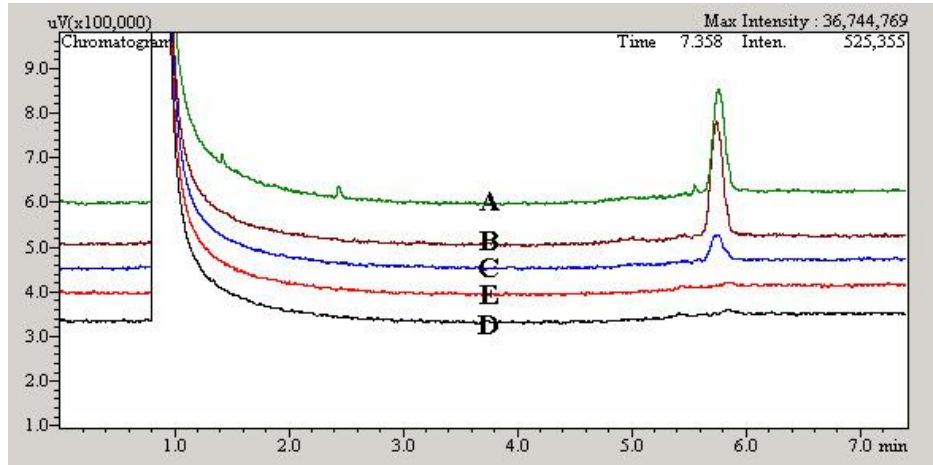


Figure 4: GC diagram of phoxim in papaya.

A phoxim 5.25 mg/kg; B sample of pericarp; C sample of the whole fruit; D and E as CK

Table1: Weight Coefficient of pesticide residues in pericarp of common pome fruit

Fruit	Quality(g)	C_p/C_T
Fuji apple	169.605±4.666	6.241
	189.425±4.748	7.963
	209.032±4.015	9.262
Pear	135.098±3.411	4.254
	146.270±3.566	4.516
	159.215±3.652	5.069
Papaya	1060.96±10.21	4.390
	1072.21±10.33	4.815
	1080.67±11.19	5.037

3.3 Formula for calculating the amount of pesticide residues in the epidermis

Although the ratio of the pesticide in pericarp and the total (C_p/C_t) ranged from 6.241 to 9.262 in apple, 4.254 to 5.069 in pears, and 4.39 to 5.037 in papayas in different groups of weight, the ratio of C_p/C_t and quality of specific group is essentially unchanged for a certain kind of fruit (Table 2). According to these results, a conversion formula of the content of pesticide residues in whole fruit and pericarp was set based on surface sampling method.

The weight of a pesticide(a) refers to the epidermis and fruit pesticide content ratio of the overall pesticide content(1).

$$a = \frac{C_p}{c_t} \quad (1)$$

a : Weight coefficient of pesticide residues in pericarp,

C_t : total pesticide content,

C_p : pesticide content in the epidermis

Coefficient A refers to the contribution of the content of epidermal pesticide to the whole fruit

$$A = \frac{a}{m} \quad (2)$$

A: surface pesticide coefficient,

a: Weight coefficient of pesticide residues in pericarp,

m: the quality of a single sample

According to the data of Weight coefficient of pesticide residues in pericarp of common pome fruits in Table1, the surface pesticide coefficient could be calculated based on formula (2). As shown in Table2, the surface pesticide coefficient A of Fuji apple, pear and papaya were 0.041, 0.031 and 0.004 respectively.

So fruit samples with pesticide residues in total mass fraction of w, in units of milligrams per kilogram (mg / kg), according to the measured pesticide residues in skin samples according to the formula (3)

$$w = \frac{w_p}{A \times m} \quad (3)$$

w - the total pesticide residues in fruits, in milligrams per kilogram (mg / kg);

w_p - the amount of pesticide residues in fruit pericarp, in milligrams per kilogram (mg / kg);

m - single sample sample quality of the sample, in grams (g);

A - surface pesticide coefficient

Table 2. Surface pesticide coefficient in common pome fruit

Fruit	A
Fuji apple	0.041
Pear	0.031
Papaya	0.004

3.4 Application of the surface sampling method on Rapid detector of Pesticides

Select the most commonly used Rapid detector of Pesticides or detector have a fast detector function ,to detect pesticides in food safety rapid detector to verify the sampling method, the results shown in Table 3. Take apples, pears and papaya samples of different quality and shape of the three groups of 15 individuals, according to this standard involve the rapid detection of pesticide residues in fruit skin sampling methods in the technical parameters, take the skin (with 2.5mm pulp) as sample, while the overall homogenized fruit sample, according to GB/T 5009.199-2003 "vegetables of organophosphorus and carbamate pesticide residues in the rapid detection," Determination of pesticide content. It can be seen from the results, obtained in accordance with the samples of fruit skin test results of pesticide residues, pesticide residues on the enzyme that the inhibition rate was significantly greater than the overall sample of inhibition of fruit (see Table 3).

Table3. Comparison of the results of rapid detectors using different sampling method

Sample	Inhibition rate (%)
Whole apple	34.6 ± 5.31
Peel of apple	84.4 ± 4.89
Whole pear	30.1 ± 5.83
Peel of pear	71.8 ± 3.26
Whole papaya	35.4 ± 2.91
Peel of papaya	91.3 ± 6.03

The result of rapid detectors using whole apple sample and surface sample demonstrated that, a surface sampling method for analysis of pesticide residues in pome fruit could increased sensitivity and reduced the average impact of whole fruit homogenates. Combination of the phoxim residues of apples, pears and papayas determined by GC-MS in different layers of the fruit and the formula for calculating the amount of pesticide residues in the epidermis amount of pesticide residues in the epidermis, an efficient surface sampling method was studied and established.

4. Discussion

According to the reseach of pesticide residues in pome fruit by layer-by-layer scan from pericarp to kernel, the previous solution of the significant differences in the contents of residues between the pericarp and entire apple was demonstrated, further more these research revealed the pesticide residues, mainly remained in the pericarp and pulp of the following 2.5mm thick, this part of the sample could be selected as sample to detect the levels of residues of pesticides, because this part of the high levels of pesticide residues, can increase rapidly screening test accuracy [18]. From the experimental results obtained, the standard design of the skin sampling methods, reduction in access to sub-samples the following 2.5mm thick skin and flesh as the test object.

Scientific and standardized sampling method can reduce the skin samples and pulp composition on the average of the impact test results, reduce false positives and improve accuracy [19]. But different types of fruit, types of pesticides and interval time after spraying should be further studied to verify the coefficient and formula. Effect of surface waxes on the persistence of pesticide in Fuji apples, pears, and papayas would effect the diffusion of pesticide to pulp.

Coefficient A refers to the epidermal skin pesticide pesticide content of the fruit content of the overall contribution of major factors affecting A include the following factors: (1) differences include different types of fruit and fruit surface area for adsorption of pesticides in different penetrate the skin of the fruit impact of pesticides [20]. (2) determination of pesticide residues in skin samples when differences in sampling methods must be standardized to ensure that the skin sample, and with the standard method for sampling the epidermis consistent with the provisions. (3) the time factor after the recommended interval of sampling.

Therefore, this surface sampling method and the corresponding formula can be applied to various methods of detection of pesticide residues, in particular, the portable rapid detection machines. As the rapid development of materials science, biology and chemistry technologies, the surface sampling method would be widely used in sensor and other equipments.

Acknowledgements

We thank the Special Fund for Agro-scientific Research in the Public Interest, the Ministry of Agriculture of the People's Republic of China (Project 201003008 and 201203046) and National Key Technology R&D Program (Project 2009BADB9B07), for financial support.

References

1. Hamers, T., P. J. van den Brink, et al. (2003). "Estrogenic and esterase-inhibiting potency in rainwater in relation to pesticide concentrations, sampling season and location." *Environ Pollut* 123(1): 47-65.
2. Fenske, R. A., A. Bradman, et al. (2005). "Lessons learned for the assessment of children's pesticide exposure: critical sampling and analytical issues for future studies." *Environ Health Perspect* 113(10): 1455-62.
3. Wang, J., L. Tuduri, et al. (2009). "Flexibility of solid-phase microextraction for passive sampling of atmospheric pesticides." *J Chromatogr A* 1216(15): 3031-7.
4. Wang, K. and J. L. Shipp (2001). "Sequential sampling plans for western flower thrips (Thysanoptera: Thripidae) on greenhouse cucumbers." *J Econ Entomol* 94(2): 579-85.
5. Ambrus, A. (2009). "Estimation of sampling uncertainty for determination of pesticide residues in plant commodities." *J Environ Sci Health B* 44(7): 627-39.
Chromatogr A 1216(15): 2972-83.
6. Yusa, V., C. Coscolla, et al. (2009). "Sampling and analysis of pesticides in ambient air." *J Chemosphere* 51(6): 509-13.
7. Neumann, M., M. Liess, et al. (2003). "A qualitative sampling method for monitoring water quality in temporary channels or point sources and its application to pesticide contamination."
8. Shen, J. and J. Liu (2009). "[Determination of carbendazim residue in orange and soil using high performance liquid chromatography]." *Se Pu* 27(3): 308-12.
9. Luan YX, Ping H, Pan LG "Residue dynamics of phoxim in pericarp, sarcocarp and kernel of apple" *IFIP AICT* 345, pp.457-464, 2011.
10. Shen, J. and J. Liu (2009). "[Determination of carbendazim residue in orange and soil using high performance liquid chromatography]." *Se Pu* 27(3): 308-12.
11. Schafer, R. B., A. Paschke, et al. (2008). "Performance of the Chemcatcher passive sampler when used to monitor 10 polar and semi-polar pesticides in 16 Central European streams, and comparison with two other sampling methods." *Water Res* 42(10-11): 2707-17.
12. Schafer, R. B., A. Paschke, et al. (2008). "Performance of the Chemcatcher passive sampler when used to monitor 10 polar and semi-polar pesticides in 16 Central European streams, and comparison with two other sampling methods." *Water Res* 42(10-11): 2707-17.
13. Schafer, R. B., V. Pettigrove, et al. "Effects of Pesticides Monitored with Three Sampling Methods in 24 Sites on Macroinvertebrates and Microorganisms." *Environ Sci Technol*.
14. Bernard, C. E., M. R. Berry, et al. (2008). "Sampling household surfaces for pesticide residues: comparison between a press sampler and solvent-moistened wipes." *Sci Total Environ* 389(2-3): 514-21.
15. Rosenheck, L., J. Cowell, et al. (2001). "Determination of a standardized sampling technique for pesticide transferable turf residues." *Bull Environ Contam Toxicol* 67(6): 780-6.

16. Teng, H., V. Kwigizile, et al. (2008). "Investigation of the AP-42 sampling method." *J Air Waste Manag Assoc* 58(11): 1422-33.
17. Li, H., E. L. Vermeirssen, et al. "Controlled field evaluation of water flow rate effects on sampling polar organic compounds using polar organic chemical integrative samplers." *Environ Toxicol Chem* 29(11): 2461-9.
18. Rabiet, M., C. Margoum, et al. "Assessing pesticide concentrations and fluxes in the stream of a small vineyard catchment--effect of sampling frequency." *Environ Pollut* 158(3): 737-48.
19. Berger, B. and C. von Holst (2001). "Pesticide residues in products of plant origin in the European Union. Sampling strategy and results from the co-ordinated EU monitoring programmes in 1996 and 1997." *Environ Sci Pollut Res Int* 8(2): 109-12.
20. Gaultier, J., A. Farenhorst, et al. (2008). "Regional assessment of herbicide sorption and degradation in two sampling years." *J Environ Qual* 37(5): 1825-36.